# STEGO

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**KÄIMISRAAM**

**WALKING FRAME**

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**EN****PREFACE**

- Thank you for the confidence shown in Kid-Man®. Every product is factory checked and leaves us in faultless condition. We wish you much success and hope that the Kid-Man® product will have a positive impact on your daily life.

**INDICATIONS**


- Kid-Man® walking frame is suitable for everyday indoor and outdoor use. All Kid-Man® walking frames are manufactured to the highest quality standards and are CE marked as medical devices. The walking frame is made of anodized aluminium tubes.

**SET UP**


- The walking frame comes fully assembled and ready to use.

**USAGE**


- Hold the padded handles and pull both sides outward until the push button locks into place. You will hear a click. When ready for use, the walking frame forms a “U” shape (Pic. No. 2). Once the frame is opened, set the height of the frame according to the user’s needs (Pic. No. 3).
- When not in use, it is recommended to fold the frame. With one hand, press the plastic holder in the middle of the frame, fold one side inward, and repeat the same procedure with the other side (Pic. No. 1). The folded frame takes up less space and is more convenient for storage.


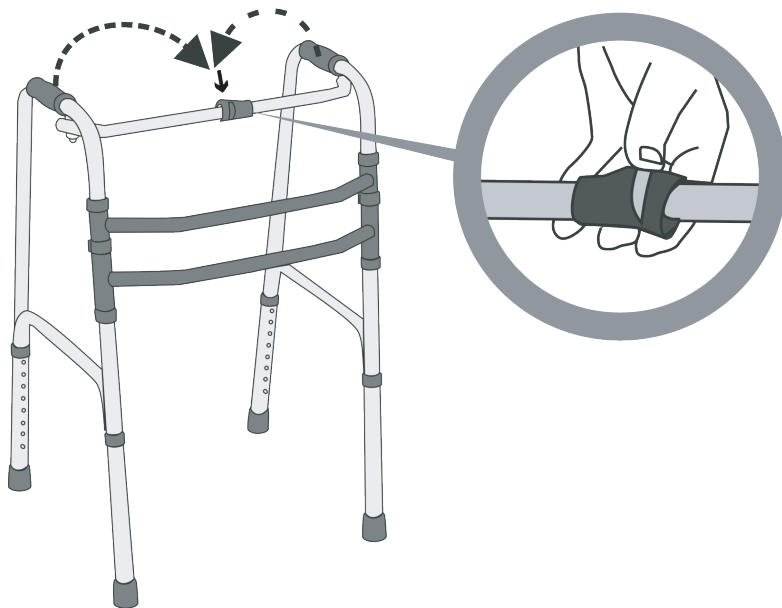
Pic. No. 1

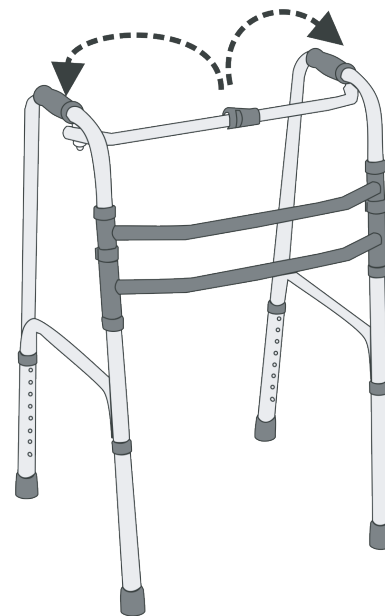
Pic. No. 2

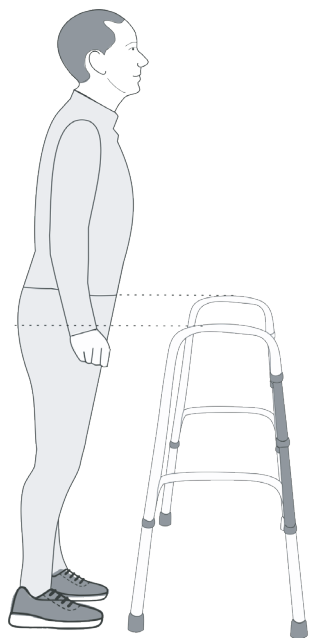Pic. No. 3

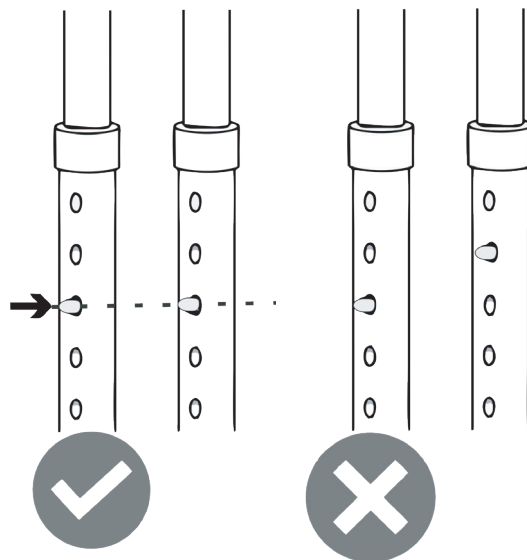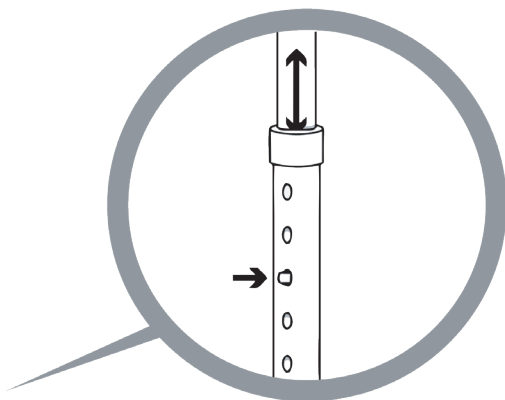• **Warning!** Before using the device, make sure all feet are set to the same height.
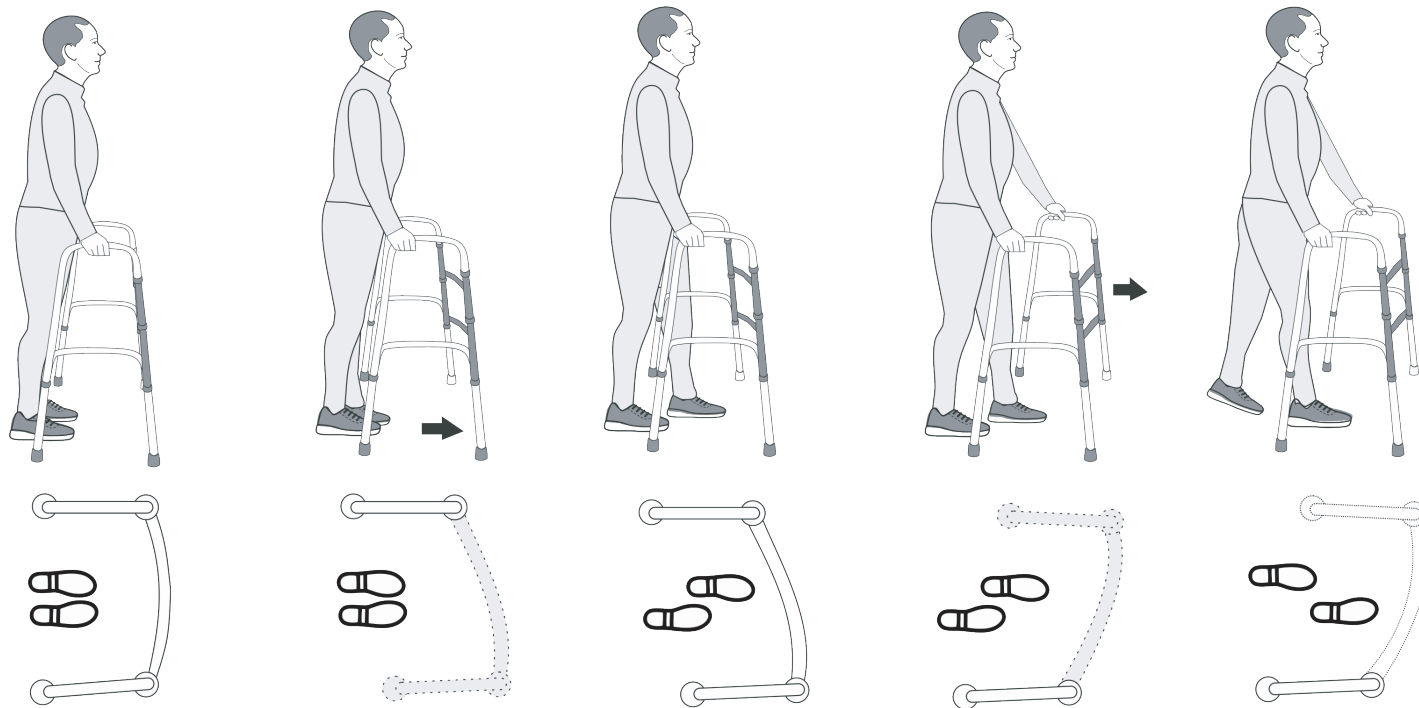
#### CARE AND MAINTENANCE

- It is recommended to clean with a damp cloth soaked in warm water. Do not use aggressive cleaning agents or solvents.

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**EN** We issue a warranty of 24 months from date of purchase. Any unauthorized modifications carried out without approval of the manufacturer will make the liability null and void. This includes exceeding the maximum user loads and using product in improper way. Warranty is only valid with fully filled warranty card. UAB „Kasko Group“ after getting claimed item from purchaser within period of 4 weeks decides whether issue is treated under warranty. If issue is not covered under warranty and customer agrees, UAB „Kasko Group“, if possible, fixes issue. All costs, including costs of transportation, should be covered by customer. Warranty repair does not extend warranty.

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